

# Contiplex<sup>®</sup> Echo

Styleted Echogenic CPNB Catheter System



Designed with echogenicity in mind.

# Contiplex® Echo

## Echogenic Catheter and Needle System for CPNB

### Styleted Threadability

Help ensure control while threading



### Flat Coil Technology

Entire catheter length constructed with flat metal coils designed to be echogenic

### Contiplex® Tuohy Ultra stimulating needle

with microdot electrode designed to provide more pronounced motor response at closer distances to nerve




Laser-etched markings along distal 2 cm, designed to be echogenic

### Patented Sideport Valve


Inject, hydrodissect and aspirate with no need to disconnect before threading catheter

### Your Choice

Available in open-tip or closed-tip configurations

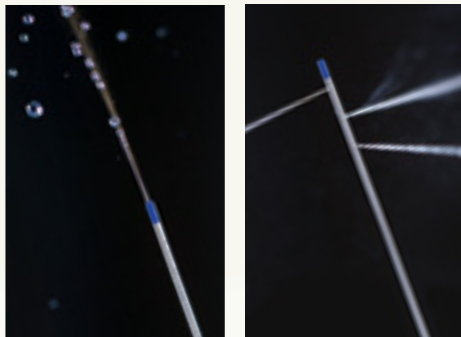

## Echogenicity in mind.

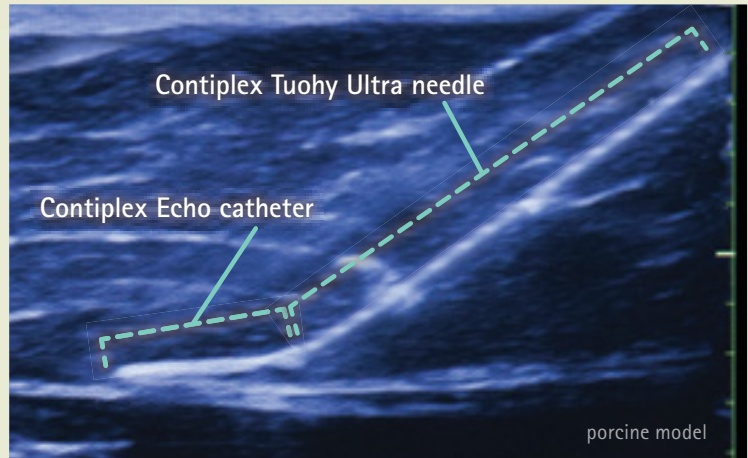

### Standard Configurations\*

Item Number	Contiplex Tuohy Ultra Needle	Contiplex Echo Catheter
331761	18G x 2"	Closed Tip
331763	18G X 4"	Closed Tip
331764	18G x 2"	Open Tip
331762	18G x 4"	Open Tip

Above sets include: Contiplex® Echo catheter, Contiplex® Tuohy Ultra needle, sideport valve and catheter connector.

### Support Tray

332103	NB400 support tray with transducer cover, ultrasound gel, tinted ChloroPrep®, ultrasound drape with 6" x 4" fenestration, fixation, and other procedural components
--------	---

\*Contiplex Echo catheters also available through B. Braun CustomPlus™ Custom Tray Program.

Schedule a demo today.

Call your B. Braun Representative at 1-800-227-2862.

[www.bbraunusa.com](http://www.bbraunusa.com)

**B | BRAUN**  
SHARING EXPERTISE