

DIALOG+[®] HEMODIALYSIS SYSTEM JUST GOT UPGRADED

With Software 9.1B Comfort Plus


+ NEW SOFTWARE 9.1B COMFORT PLUS FEATURES

- ▶ **Automatic Blood Pressure Module: Designed for Comfort**
 - » Ability to reset without turning machine on/off
 - » Measures during inflation (only with M2015 ABPM)
- ▶ **Increased Regulation:** additional sensor added to monitor dialysis fluid temperature
- ▶ **Optimized Flow Control:** added balance chamber calibration
- ▶ **Optimized Dialysate Flow:** air is purged seamlessly from DF filter every 150 seconds
 - » Alarms if air separator sensor opens more than 2 times in one hour

BP
COMFORT

+ PLUS, TRUSTED FEATURES YOU CAN RELY ON

- ▶ **Adimea™ Dose Monitor:** Continuous real-time measurement of dialysis dose
 - » Improvements to the prediction curve with 9.1B Comfort Plus software
 - » Change parameters to meet Kt/V goals during treatment
 - » Only input required is patient's weight and target Kt/V

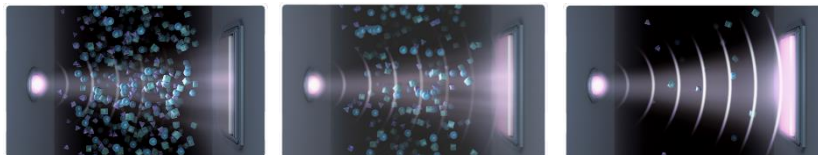

- ▶ **Automatic Programming**
 - » Program to automatically turn on/off
 - » Program to run specific selections, including disinfection and start preparation
- ▶ **Event Timer:** a reminder for individual events
- ▶ **Patient Card Reader**
- ▶ **Profiling Capability:** includes 23 pre-programmed UF profiles, such as Na, bicarb, dialysate flow, temperature, and heparin
- ▶ **Rotating 15" LCD Touchscreen Monitor:** user friendly icons
- ▶ **Stand-by Mode:** prevents waste by stopping dialysate flow at the end of preparation
- ▶ **Technical Services Mode:** test internal components and set system parameters through monitor



+ PLUS, RELIABLE CUSTOMER SUPPORT

- ▶ Free clinical webinar training
- ▶ Technical training classes
- ▶ 24/7 technical phone support
- ▶ Premium service packages