

Contiplex® Stim Stimulating Catheter Nerve Block System

Simplifies and helps increase procedural
efficiency and safety



Enhanced catheter
with added stiffness
and increased flow rate

Optimized for ease of use—change stimulating current from
needle to catheter with the flick of a switch.

B. Efficient.

B | BRAUN
SHARING EXPERTISE

The Contiplex® Stim System with 55cm Catheter

Simplifies and helps increase procedural efficiency and safety

Helps make catheter placement as easy as possible

The system consists of the Contiplex Tuohy Needle, 55cm catheter, Stimuplex® HNS12 Stimulator and Switch.

- Sets up easily before the procedure to help you keep focus on the patient during the procedure
- System can be used with or without the switch*



*Switch can only be used with Stimuplex HNS11 and HNS12 Nerve Stimulators

Reduce procedural steps



Plug needle and catheter into switch box before procedure




Change stimulation current from needle to catheter with the simple flick of a switch

B. Safe.

Enhance safety and efficiency with needle and patented sideport valve



Pre-load catheter in sideport valve before procedure to reduce chance of catheter leaving sterile field

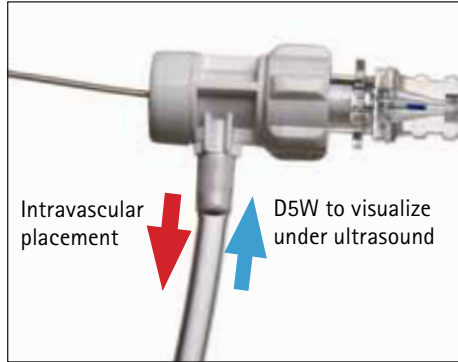

Aspirate and inject with **pre-attached** sideport for convenience




Microdot Electrode tip on needle focuses current more precisely for optimal motor response


55cm catheter provides easy handling and is clearly visible under ultrasound*



Stylet reinforces catheter for added stiffness and easier threading

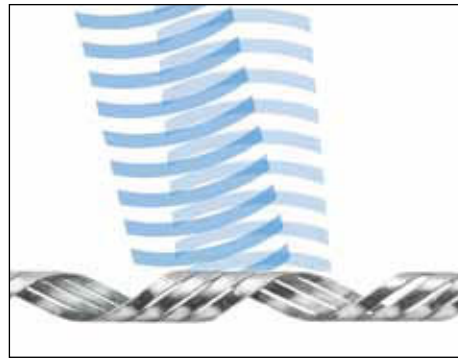

Four flat-wire coils maximize reflective surface to provide clear visibility under ultrasound

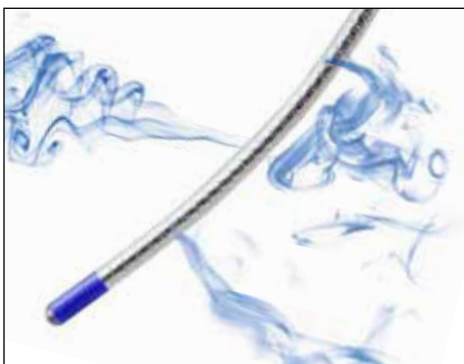

Increased flow rate from shorter catheter



Fully enclosed, flexible tip helps reduce trauma and enhances glide through needle

* Stuart Grant, MD, Duke University

Optimize your peripheral nerve block procedures with these innovative products:

Stimuplex® HNS12 Nerve Stimulator

- Provides increased accuracy through constant visual feedback regarding quality of motor twitches being elicited.
- Received highest score when evaluated by anesthesiologists with 12 other stimulators.¹
- Easy operation through intuitive menus and quick access keys.



Contiplex® Stim Nerve Block Tray

With **NEW** tinted ChloroPrep® and 55cm catheter in polybag



Sharing Expertise®

It's more than a philosophy.
It's a promise.

Sharing Expertise includes listening to the wisdom of colleagues, customers and clinicians. It's more than a philosophy. It's our promise to the healthcare community and ourselves to continuously seek knowledge.

To expand it. To apply it. To share it.
All to help improve and save lives.

This is Sharing Expertise.

This is B. Braun.

Schedule a demo today.

Call your B. Braun Representative at **1-800-227-2862** or visit us at **www.contiplexstim.com**

Ordering Information

REF #	Product Code	Description	Qty./case
331755	SCNB2TK	Contiplex Stim Nerve Block Tray (2 in. needle): Contains 2" Contiplex Tuohy needle, Contiplex Stim catheter in dispenser, filter, catheter connector, PNB label and syringes, drape, drugs, needles, tinted ChloroPrep.	10
331756	SCNB4TK	Contiplex Stim Nerve Block Tray (4 in. needle): Contains 4" Contiplex Tuohy needle, Contiplex Stim catheter in dispenser, filter, catheter connector, PNB label and syringes, drape, drugs, needles, tinted ChloroPrep.	10
331757	SCNB2	Contiplex Stim Nerve Block Set (2 in. needle): Contains catheter, needle, sideport, PNB label, catheter connector, filter	10
331758	SCNB4	Contiplex Stim Nerve Block Set (4 in. needle): Contains catheter, needle, sideport, PNB label, catheter connector, filter	10
333543	STIMSWITCH	Stimuplex Switch: Nerve stimulator current selector between needle/catheter — Only for use with Stimuplex HNS11 and HNS12.	1

B | BRAUN
SHARING EXPERTISE

www.bbraunusa.com

1) Jochum, D., et al. An objective assessment of nerve stimulators used for peripheral nerve blockade. Anesthesia 2006, 61, pages 557-564
Rx only. © 2009 B. Braun Medical Inc., Bethlehem, PA. All rights reserved. ChloroPrep is a registered trademark of Cardinal Health, Inc.